EN DeVilbiss® SmartLink® Module Instruction Guide, Model DV5M

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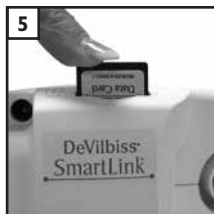
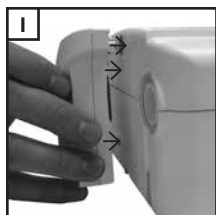
INTRODUCTION

The DeVilbiss SmartLink Module is an optional accessory to the DV5X series CPAP devices. This SmartLink device was supplied to you by your homecare provider as part of your overall positive airway pressure therapy for obstructive sleep apnea (OSA). SmartLink is a helpful tool for your physician and/or homecare provider because it allows them to see the effectiveness of the therapy and how frequently and consistently you use the device.

The information SmartLink collects helps bring to light potential problems or challenges you are facing that your clinician might be able to help you with in getting used to therapy.

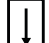
SmartLink is an advanced technology module that works together with your CPAP (Continuous Positive Airway Pressure) device.

This system is unique in healthcare because many times a clinician will prescribe a therapy and send you home with it and never really know how effective it is. SmartLink is a great benefit to you because it helps your clinician keep track of what's going on with your condition without you having to make extra trips to the clinician's office.

Indications For Use

The DeVilbiss SmartLink System can only be used in conjunction with the DV51, DV53, DV54, DV55, DV56, and DV57 Series CPAP Systems for follow up of obstructive sleep apnea in patients weighing above 30 kg on nasal CPAP therapy.

SYMBOL DEFINITIONS

	Attention - Consult Instruction Guide		Data Port Input/Output
	Type BF equipment- applied part		Memory Card Insertion Orientation
	Class II electrical protection- double insulated		Oximeter Connection Port
	This device contains electrical and/or electronic equipment that must be recycled per EC Directive 2002/96/EC- Waste Electrical and Electronic Equipment (WEEE)		

IMPORTANT SAFEGUARDS

Please read all instructions carefully before use. Important information is highlighted by the following terms:

DANGER Urgent safety information for hazards that will cause serious injury or death.

WARNING Important safety information for hazards that might cause serious injury.

CAUTION Information for preventing damage to the product.

NOTE Information to which you should pay special attention.

DANGER

- Electric Shock Hazard-Do not use while bathing.
- Electric Shock Hazard-Do not immerse this device into water or any other liquid.
- Electric Shock Hazard-Do not attempt to open or remove the enclosure. There are no user-serviceable components inside. Contact your home care provider if your device needs servicing. Opening or tampering with the product will void the warranty.

WARNING

- Use the SmartLink Module only for its intended use as described in this manual.
- Be sure to read and understand all safety instructions supplied with your SmartLink Module device.
- Use only accessories recommended by DeVilbiss.
- Use only data cards provided by your homecare provider. Other Secure Digital (SD) cards may not be compatible with the system or may not perform properly.
- Do not expose the device to water.
- If water has been spilled onto the SmartLink Module or it has been submerged into water, unplug the power cord from the power source. Allow the device to completely dry before use.
- If the device has been dropped, refer to the troubleshooting guide for instructions. Contact your home care provider for evaluation of the device.
- To avoid electric shock, always unplug power cord from wall outlet power source when performing cleaning.
- Do not insert any foreign objects into any openings.

CAUTION—The circular data port and oximeter connectors located on the side of the SmartLink Module are used to attach accessories to the device. The connectors must only be used with accessories approved for use by DeVilbiss. Do not attempt to attach any other device to these connectors as it may damage the CPAP, SmartLink Module, and/or the device you are attempting to attach.

Be sure that you read and understand the information in this instruction guide. The information in this guide will help you understand the operation of the Smartlink Module and how to properly maintain it. By following these instructions and the advice of your physician, your SmartLink Module will become an effective addition to your CPAP therapy.

SAVE THESE INSTRUCTIONS

SET UP

Important Parts Of Your SmartLink Module

1. Data card slot
2. Release button
3. Secondary accessory data connection port
4. Oximeter connection port
5. Primary data connector
6. Locking tabs (3 tabs)
7. Instruction Guide (not shown)

Connection To The CPAP

NOTE—Your CPAP may already have been provided to you with the SmartLink Module attached. If you have received your module separately from the CPAP, follow the steps below to properly attach the module to the CPAP.



1. Locate the primary data connector on the back of the SmartLink Module. Align and engage the primary data connector into the circular DIN connector on the back of the CPAP. At the same time, align and engage the three locking tabs on the module with the three locking recesses on the back surface of the CPAP **Figure 1**.
2. **Firmly** press the module onto the back of the CPAP until the module is fully secured and locked into position **Figure 2**. An audible click may be heard confirming that the module is properly attached to the CPAP. Visually confirm that there is not a gap between the module and the CPAP.

Removal From The CPAP

To remove the module from the CPAP:

1. Fully depress the release button on the side of the module.
2. Pull the module straight back and away from the CPAP.

CAUTION—Do not pull up on the module as this may damage the data connector.

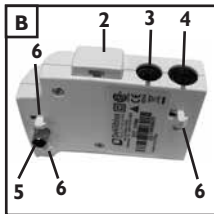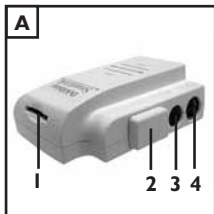Optional Oximeter Connection

The SmartLink Module provides oximeter connection capability to allow an approved Nonin XPOD oximeter to be used with the CPAP/module system. Connect the oximeter cable into the oximeter connection port located on the side of the module. If you are using an oximeter, please refer to the assembly and operation instructions that came with the oximeter **Figure 3**.

You can see your oximeter readings by using the SmartLink menu (see SmartLink Menu Structure) at any time.

Message shown when the oximeter is first connected and logging starts. This message will remain on display for five seconds. The pulse oximeter menu (SpO₂/PR) will then be displayed.

Oximeter
Connected

Message shown when the oximeter is removed. This message will remain on display for five seconds. Multiple notifications will occur until logging stops.

Oximeter
Removed

SmartLink Menu Structure

To display or remove SmartLink menu, press Left and Right Arrow Keys together. Use Arrows to navigate through the menu structure as shown (will follow dotted line).

Navigation:

- ▶ Right Arrow - navigates right through menu items.
- ◀ Left Arrow - navigates left through menu items.
- ▲ Up Arrow - navigates up through menu items.
- ▼ Down Arrow - navigates down through menu items.

To Exit SmartLink menu (choose one):

1. Turn blower on or off
2. Press Left and Right Arrow Keys together.

SmartLink Menu Press ◀▶ To Exit	Volume 100% (Adjust ▲▼)	Serial Number CM123456	Firmware Version V0.06 25/04/2008	SP02 100 PR: 101
------------------------------------	----------------------------	---------------------------	--------------------------------------	---------------------

OPERATING INSTRUCTIONS

Navigating Menus and Onboard Reporting

While the SmartLink Module is attached, it will collect important data about how the device is being used and what kind of therapy it is delivering for you. Much of the information that is collected is available to you through the onboard menu on your CPAP display. The previous table shows how to navigate through the menus.

When the SmartLink is connected correctly, the display on the CPAP will show this message. **NOTE**—if a data card is inserted in the module before attaching the module to the CPAP, then card detection messages may be displayed.

SmartLink Module
Ready for Use

To enter the SmartLink menu, press the LEFT and RIGHT Arrow keys on the CPAP keypad at the same time. To exit the SmartLink menu, press the same key sequence again or simply turn the CPAP off and on.

SmartLink Menu
Press ◀▶ To Exit

Press the UP or DOWN arrow key on the CPAP keypad to adjust the volume for the audible alert notification.

Volume
100% (Adjust ▲▼)

This menu displays the serial number of the SmartLink module connected to the CPAP.

Serial Number
CM123456

This menu displays the operating firmware version of the SmartLink module connected to the CPAP.

Firmware Version
V0.06 25/04/2008

This menu displays the pulse oximetry (SpO2) and pulse rate (PR) information when a pulse oximeter is used with the system. If a pulse oximeter is not connected to the module "****" will be displayed for the SpO2 and PR values.

SpO2: 100
PR: 101

Data Retrieval

Your clinician may have a standard protocol to follow-up on your therapy to find out how the device is working and to identify any problems or issues you may be having. Your home care provider will mail a SmartLink Data Card to you at certain intervals and ask you to download the data onto the card and mail it back. It is a very simple process:

1. Ensure that the SmartLink module is properly connected to the CPAP and that the CPAP is connected to an appropriate power source.
2. Orient that card as shown in **Figure 4**. Insert and press the card into the card slot located on the top of the module until the card clicks into place **Figure 5**.
3. Ensure that the card is detected (view message or audible alert) and allow the data to transfer to the card. Depending on the amount of data, this could take several minutes. Do not remove the card from the module until complete. Observe the message(s) on the display to confirm that the transfer was successful (see Messages below). If necessary, contact your provider and report any errors that may have occurred.
4. If instructed by your provider, eject the card from the card slot by pressing down on the card until it is released and remove the card from the slot.
5. Place the card in the provided mailer and return it to the provider.

Messages

The following messages could be displayed on the CPAP during the data retrieval process:

Message shown when a card is first inserted into the SmartLink Module and data is being transferred. This message will remain on the display until transfer is completed.

Card Detected
Please Wait . . .

Message shown when a card is first inserted into SmartLink Module and there is no new data to transfer. Will remain on display for ten seconds.

Card Detected
Ready For Use

Message shown after completion of card data transfer. This message will remain on the display for ten seconds.

Card Transfer
Successful

Message shown if an error is encountered during card data transfer. This message will remain on the display until acknowledged by pressing any key on the keypad or a 15 minute timeout occurs.

Card Error
Contact Provider

Message shown if an error is encountered while updating device settings. This message will remain on the display until acknowledged by pressing any key on the keypad or a 15 minute timeout occurs. **NOTE**—If this message is displayed, contact your clinician.

Setting Error
Contact Provider

Message shown after removing the card from the module. This message will remain on the display for ten seconds.

Card Removed
Send to Provider

Message shown when a new module with mismatched data is detected.

This message will remain on the display until the module is removed.

Contact your provider to notify them of the situation and for instructions on how to resolve the issue. If this situation is not resolved, the module will not be able to record usage and therapy data. However, although this message will remain on the display with the module attached, the CPAP will still be functional and provide therapy.

Setup Required
Contact Provider

Other Messages

In the unlikely event that a fault should occur, this message will be shown when a module fault condition is encountered (Clock Error, Memory Error, Setting Error, etc).

E00 Module Fault
Contact Provider

Message shown when firmware is being updated. This message will remain on the display until the update is complete. Occasionally a firmware update may be provided for your device via SmartLink Data Card. The firmware update process is automatic and may take a few seconds to complete. Once the update is complete, the system will return to normal operation.

Firmware Update
Please Wait . . .

ACCESSORIES

Oximeter Kit.....	DV5M-697
Oximeter Finger Clip Sensor	8000AA
Oximeter Ear Clip Sensor	8000Q
Oximeter Finger Soft Sensor, Medium	8000SM
Oximeter Finger Soft Sensor, Large.....	8000SL

CLEANING INSTRUCTIONS

WARNING

To avoid electric shock, always unplug power cord from the wall outlet power source.

CAUTION-Never rinse or place the device in water. Never allow liquids to get into or around any of the ports or openings; doing so will result in product damage. If this occurs, do not use the device. Contact your provider for service.

NOTE-The module should be cleaned as needed.

1. Turn off the CPAP and unplug it from the power source.
2. Use a soft damp cloth to clean the surfaces of the device.
3. Allow the device to dry completely before plugging it into a power source.

SERVICE INSTRUCTIONS

The SmartLink Module does not require regular servicing.

TROUBLESHOOTING

Problem	Probable Cause	Remedy
Module is not recognized when connected to the CPAP.	1. The CPAP is not connected to a power source.	1. Connect the CPAP to an appropriate power source and verify that the module is recognized.
	2. Module is not fully attached to the CPAP.	2. Remove the module from the CPAP, realign the primary data connector with the circular DIN connector of CPAP, and press firmly until the module is locked into position. Be sure there is not a gap between the module and the CPAP.
Card not recognized when inserted into the module card slot.	1. Card not oriented properly.	1. Refer to the symbol on the top of the module or Figure 4 in this guide for proper card orientation.
	2. Card not fully inserted in the module card slot.	2. Press the card fully into the card slot until the card clicks into position.
	3. Module is not fully attached to the CPAP.	3. Remove the module from the CPAP, realign the primary data connector with the circular DIN connector of CPAP, and press firmly until the module is locked into position. Be sure there is not a gap between the module and the CPAP.
Card error.	1. Wrong card type.	1. Only use a card supplied by your provider. The card must be an SD card with FAT16 format.
	2. Card removed during data transfer.	2. Card must not be removed during data transfer. Data transfer interruption may corrupt the card and prevent the card from working correctly. Reinsert card and allow data transfer to continue until completed.
	3. Card corrupt or faulty.	3. Return card to provider for replacement.
Settings error.	1. Card removed during data transfer.	1. Card must not be removed during data transfer. Data transfer interruption may corrupt the card and prevent the card from working correctly. Reinsert card and allow data transfer to continue until completed. If this error is not corrected, contact provider for resolution.
	2. Card corrupt or faulty.	2. Return card to provider for replacement. If this error is not corrected, contact provider for resolution.
Setup required.	1. Module not set up properly.	1. Notify provider of the situation and follow provider instructions.
Oximeter data (SpO ₂ and PR) not displayed.	1. Oximeter not properly connected to the SmartLink module.	1. Verify connection between the oximeter connector and module oximeter connection port.
	2. Oximeter probe not properly attached to user.	2. Reference the operating instructions that came with oximeter.

Problem	Probable Cause	Remedy
CPAP is not fully functional.	1. Module may be defective.	1. Remove the module from the CPAP and verify that CPAP is fully functional. If it functions properly without the module connected, the module may be defective. Contact provider. If it still is not fully functional, consult the CPAP instruction guide for troubleshooting.

DEVILBISS SMARTLINK SPECIFICATIONS

Size: 3.9" H x 3.1" W x 1.4" D (9.9cm H x 7.9cm W x 3.6cm D)

Weight: 0.30 lbs. (0.14 Kg)

Electrical Rating:

Maximum Power Consumption (SmartLink Module only): 4 W

Voltage and Current: Power supplied by the DV5X Series Device

Class II Equipment; Type BF Applied Parts; Continuous Operation

Operating Conditions:

Temperature Range: 41°F to 104°F (5°C to 40°C)

Humidity Range: 0 to 95% R.H. non-condensing

Atmospheric Conditions: Sea level to 9,000' (2,743 M)

Transportation and Storage Conditions:

Temperature Range: -40°F to 158°F (-40°C to +70°C)

Humidity Range: 0 to 95% R.H. non-condensing

DEVILBISS GUIDANCE AND MANUFACTURER'S DECLARATION

WARNING

Medical Electrical Equipment needs special precautions regarding EMC and needs to be installed and put into service according to the Electromagnetic Compatibility [EMC] information provided in the accompanying documents.

Portable and Mobile RF Communications Equipment can affect Medical Electrical Equipment.

The equipment or system should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the equipment or system should be observed to verify normal operation in the configuration in which it will be used.

NOTE– The EMC tables and other guidelines provide information to the customer or user that is essential in determining the suitability of the Equipment or System for the Electromagnetic Environment of use, and in managing the Electromagnetic Environment of use to permit the Equipment or System to perform its intended use without disturbing other Equipment and Systems or non-medical electrical equipment.

Guidance and Manufacturer's Declaration – Emissions All Equipment and Systems

This device is intended for use in the electromagnetic environment specified below. The customer or user of this device should assure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic Enforcement – Guidance	
RF Emissions CISPR 11	Group 1	This device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.	
RF Emissions CISPR 11	Class B Radiated and Conducted Emissions	This device is suitable for use in all establishments including domestic, and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.	
Harmonics IEC 61000-3-2	Class A		
Flicker IEC 61000-3-3	Complies		
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	±6kV contact ±8kV air	±6kV contact ±8kV air	Floors should be wood, concrete or ceramic tile. If floors are synthetic, the relative humidity should be at least 30%.
Electrical Fast Transient/burst IEC 61000-4-4	±2kV on AC Mains	±2kV on AC Mains	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1kV Differential ±2kV Common	±1kV Differential ±2kV Common	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	>95% Dip for 0.5 Cycle 60% Dip for 5 Cycles 30% Dip for 25 Cycles >95% Dip for 5 Seconds	>95% Dip for 0.5 Cycle 60% Dip for 5 Cycles 30% Dip for 25 Cycles >95% Dip for 5 Seconds	Mains power quality should be that of a typical commercial or hospital environment. If the user of this device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or battery.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Power Frequency 50/60Hz Magnetic Field IEC 61000-4-8	3A/m	3A/m	Power frequency magnetic fields should be that of a typical location in a typical commercial or hospital environment.
Conducted RF IEC 61000-4-6	3 Vrms from 150 kHz to 80 MHz	$V1 = 3 \text{ Vrms}$	Portable and mobile RF communications equipment should be separated from the device by no less than the recommended separation distances calculated/listed below: $D = (3.5/V1)\sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	$E1 = 3 \text{ V/m}$	$D = (3.5/E1)\sqrt{P}$ 80 to 800 MHz $D = (7/E1)\sqrt{P}$ 800 MHz to 2.5 GHz Where P is the maximum power rating in watts and D is the recommended separation distance in meters. Field strengths from fixed transmitters, as determined by an electromagnetic site survey, should be less than the compliance levels (V1 and E1). Interference may occur in the vicinity of equipment containing a transmitter.
<p>For transmitters rated at a maximum output power not listed above, the recommended separation distance D in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.</p> <p>Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.</p> <p>Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.</p>			

Recommended Separation Distances Between Portable and Mobile RF Communications Equipment and this device. This device and system are NOT Life-Supporting

This device is intended for use in the electromagnetic environment in which radiated disturbances are controlled. The customer or user of this device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF Communications Equipment and the device as recommended below, according to the maximum output power of the communications equipment.

Maximum Output Power (Watts)	Recommended Separation Distances for the device (meters)		
	150 kHz to 80 MHz $D=(1.1667)\sqrt{P}$	80 to 800MHz $D=(1.1667)\sqrt{P}$	800 MHz to 2.5 GHz $D=(2.3333)\sqrt{P}$
0.01	0.11667	0.11667	0.23333
0.1	0.36894	0.36894	0.73785
1	1.1667	1.1667	2.3333
10	3.6894	3.6894	7.3785
100	11.667	11.667	23.333

For transmitters rated at a maximum output power not listed above, the recommended separation distance D in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

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DeVilbiss Healthcare Pty. Limited

15 Carrington Road, Unit 5
Castle Hill NSW 2154
AUSTRALIA
+61-2-9899-3144

DeVilbiss Healthcare SAS

13/17, Rue Joseph Priestley
37100 Tours
FRANCE
+33 (0) 2 47 42 99 42

DeVilbiss Healthcare GmbH

Kamenzer Straße 3
68309 Mannheim
GERMANY
+49-621-178-98-230

A-DV5M-2 Rev. E

DeVilbiss Healthcare LLC • 100 DeVilbiss Drive • Somerset, PA 15501 • USA
800-338-1988 • 814-443-4881 • www.DeVilbissHealthcare.com

DeVilbiss® and SmartLink® are registered trademarks of DeVilbiss Healthcare.

© 2013 DeVilbiss Healthcare LLC. 02.13 All Rights Reserved.

A-DV5M-2 Rev E

EN ENGLISH

Congratulations on your purchase of a SmartLink II Module with Bluetooth (DV5MB). Special information relating to the DV5MB is below; please refer to the included SmartLink Module (DV5M) instruction guide for all other information including use instructions.

FCC and Industry Canada

Host manufacturer: DeVilbiss Healthcare LLC

Limited Host Brand name: SmartLink II

Host model number: DV5MB

FCC ID: 2AA2G-001

IC: 3941A-001

This device complies with part 15 of the FCC Rules and with Industry Canada license-exempt RSS standard(s).

Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation of the device.

The output power of this device is less than 20mW. The SAR test is not required. When using it, ensure that the antenna of the device is at least 20cm away from all persons.

This Class B digital apparatus complies with Canadian ICES-003.

NOTE: *The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications or changes to this equipment. Such modifications or changes could void the user's authority to operate the equipment.*

NOTE: *This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:*

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CE Declaration

Hereby, DeVilbiss Healthcare LLC, declares that this SmartLink II Compliance Module with Bluetooth is in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC. The declaration of conformity may be consulted at www.devilbisshealthcare.com/files/CE_Declaration.pdf.

The Bluetooth[®] word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by DeVilbiss is under license. Other trademarks and trade names are those of their respective owners.

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